



Med Sled First Call

SIMPLIFY YOUR REMOVAL. Is your biggest challenge first call removals? Do you have to deal with bariatric deceased, small hallways and doors, or stairwells? The Med Sled First Calls helps you and your staff with the most challenging of first calls removals.



Standard

The Med Sled First Call Sled is simple and easy to deploy. The standard 90cm sled is designed for the average size deceased and allows anyone to transport a deceased person up to 3 times their own weight.


Size rolled	23cm X 92cm
Size open	92cm x 220cm
Weight	4.5kg
Maximum load	450kg



Bariatric

The Med Sled First Call Bariatric 120cm wide sled is 30cm wider than our standard 90cm Sled to accommodate larger girth of the deceased up to 450kg.

Size rolled	120cm X 25cm
Size open	120cm x 220cm
Weight	5.5kg
Maximum load	450kg

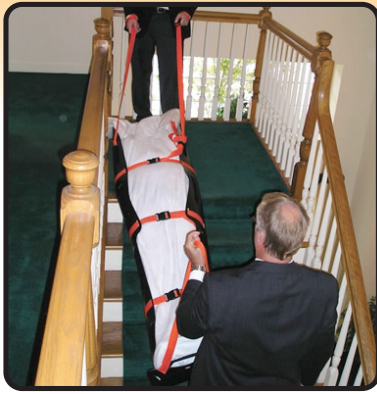

KEY BENEFITS

- Simple and easy to use
- Better than a Backboard
- Stairwell braking system reduces your risk of back injury
- Secures limp weight
- Ideal for tight spaces and multi-story homes
- Standard and Bariatric versions
- Decon Capable
- Cot Compatible
- Holds up to 450kg.
- Cost effective

Features

Leave Your Cot At The Door

The Med Sled First Call is lightweight and easy to carry to difficult situations. Stored in a 23cm X 90cm (120cm for Bariatric Sled) shoulder carry bag and weighing just 4.5kgs, throw the sled over your shoulder and take it where you need it most.

Loading The Body

Even in the most difficult situations, the First Call is designed to enable just two people to load even your heaviest removals. The sled can help protect against lower back strain.

Securing The Body

Managing and securing awkward and difficult body types is made easy with the Med Sled First Call. The sled will cocoon up and around the deceased, to provide a secure means of transport and eliminate the risk of weight shifting during the removal.

Transporting The Body

The Med Sled First Call will make transporting the body safer and easier regardless of the size of the home. The sled can be carried out like a stretcher if space permits, but in extreme cases involving large bodies and tight spaces such as trailers, slide the sled through narrow hallways and around tight corners. You can provide a dignified removal, whilst reducing risk, injury to your staff and damage to the home.

Going Down Stairs

Carrying a body down a flight of steps is a difficult and sometimes dangerous process. The end of the sled is curved to allow a smooth descent, and in more difficult removals, deploy our stairwell braking system to ensure a safe and controlled process.

Loading Onto Your Cot

The sled will cocoon the body, so once tightened, it will only be as large as the deceased and will easily fit onto your standard cot.

Easy To Clean

There are no holes or corners in the sled, so simply use your standard disinfectant to spray and wipe away any fluids.